**IMPORTANT:** Before first use, remove label from the display.

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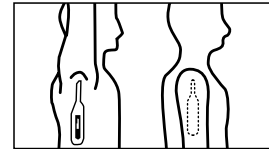
Fig. 8

## 1 YEAR LIMITED WARRANTY (Continued)

**C.** This warranty does not cover damage resulting from any unauthorized attempts to repair or from any use not in accordance with the instruction manual. Call us toll-free at 1-800-477-0457 or visit: [VicksThermometers.com/contact-us](http://VicksThermometers.com/contact-us). Please be sure to specify a model number.

**NOTE: IF YOU EXPERIENCE A PROBLEM, PLEASE CONTACT CONSUMER RELATIONS FIRST OR SEE YOUR WARRANTY. DO NOT RETURN THE PRODUCT TO THE ORIGINAL PLACE OF PURCHASE. DO NOT ATTEMPT TO OPEN THE DEVICE HOUSING YOURSELF, DOING SO MAY VOID YOUR WARRANTY AND CAUSE DAMAGE TO THE PRODUCT OR PERSONAL INJURY.**

## OPERATING YOUR VICKS RAPIDREAD™ DIGITAL THERMOMETER

**Tip:** For fastest measurements, complete steps 1 and 2 before inserting the thermometer.

### 1. TURNING THE UNIT ON

- Turn unit on by pressing the On/Off button (Fig. 1).
- A beep will sound and the display will flash Red/Yellow/Green for 1 second.
- Wait for unit to get ready
  - The display will perform a self-test (Fig. 2).
  - The display will show last measurement for 2 seconds (Fig. 3).
  - After last measurement is displayed, the display will enter ready mode (Fig. 4).

### 2. SELECTING A TEMPERATURE MEASUREMENT SITE

- After the last memory is displayed, the thermometer enters the site selection mode.
- The site selection icons will cycle from oral to underarm to rectal. (Fig. 4).
- When the icon displays the desired measurement site, press and release the blue SITE SELECTION button (Fig. 5).
- The chosen site will flash for 3 seconds and then display the ready screen (Fig. 6). You may now take a reading at the site you selected.
- You may choose an alternate site by pressing the site selection button again (Fig. 5) at any time during the 3 second confirmation period.

**WARNING:** For accurate readings, you must match the site selected on the thermometer with the measurement site on the patient.

Site	Temperature Ranges and Meaning		
Oral	96.4 - 99.4°F	99.5 - 101.3 °F	101.4+ °F
	35.8 - 37.4 °C	37.5 - 38.5 °C	38.6+ °C
Under-arm	96.2 - 99.2 °F	99.3 - 101.1 °F	101.2+ °F
	35.7 - 37.3 °C	37.4 - 38.4 °C	38.5+ °C
Rectal	97.3 - 100.3 °F	100.4 - 102.2 °F	102.3+ °F
	36.3 - 37.9 °C	38.0 - 39.0 °C	39.1+ °C
	Normal (Green)	Fever (Yellow)	High Fever (Red)

### 3. TAKING A TEMPERATURE MEASUREMENT

#### Oral Use:

- YOU MUST COVER THE THERMOMETER'S TIP WITH A SINGLE-USE, DISPOSABLE PROBE COVER.** Probe cover not included. Single-use, disposable probe covers for digital stick thermometers are available for purchase at retail chain stores and pharmacies.
- Place thermometer under the tongue. It is important to place the tip of the thermometer well under the tongue to get a good, accurate reading. (Fig. 7) Place sensor tip in areas marked with "✓". Keep mouth closed and sit still to help ensure an accurate measurement.
- When the measurement is complete, two beeps will sound. You may then remove the thermometer. The thermometer will then beep with the displayed reading.

d. Dispose of used probe cover (not included) and wash the tip as recommended.

**Under-the-arm (Axillary) Use:** This method is used for babies and young children.

- YOU MUST COVER THE THERMOMETER'S TIP WITH A SINGLE-USE, DISPOSABLE PROBE COVER.** Probe cover not included. Single-use, disposable probe covers for digital stick thermometers are available for purchase at retail chain stores and pharmacies.
- Dry the underarm with a towel before measuring.
- Place the probe tip under the arm so the tip is touching the skin and position your arm next to your body (Fig. 8) . With a young child, it is sometimes helpful to hug the child to keep their arm next to their body. This ensures that the room air does not affect the reading.
- When the measurement is complete, two beeps will sound. You may then remove the thermometer. The thermometer will then beep with the displayed reading.

**Note:** This method may produce a temperature 1°F lower than a temperature taken orally.

d. Dispose of used probe cover (not included) and wash the tip as recommended.

**Rectal Use:** Commonly used for babies, young children, or when it is difficult to take an oral or underarm temperature.

- YOU MUST COVER THE THERMOMETER'S TIP WITH A SINGLE-USE, DISPOSABLE PROBE COVER AND LUBRICATE WITH A WATER SOLUBLE JELLY FOR EASIER INSERTION.** DO NOT use petroleum jelly. Read the label of the lubricant to be sure you are using the right kind. Probe cover not included. Single-use, disposable probe covers for digital stick thermometers are available for purchase at retail chain stores and pharmacies.
  - Lay the individual being measured on his/her side. If the individual is an infant, the proper position for a baby is to lay on its stomach with legs hanging down, either across your knees or at the edge of a bed or changing table. This positions the infant's rectum for safe and easy insertion of the thermometer.
  - Gently insert the tip of the thermometer NO MORE THAN 1/2 inch into the rectum. If you detect resistance, stop. Hold the thermometer in place during measurement.
  - When the measurement is complete, two beeps will sound. You may then remove the thermometer. The thermometer will then beep with the displayed reading.
- Note:** This method produces a temperature approximately 1°F higher than a temperature taken orally.
- Dispose of used probe cover (not included) and wash and disinfect the tip as recommended.

### 4. SET A PREFERRED TEMPERATURE MEASUREMENT SITE

To save a preferred temperature measurement site, so that the preferred site will be automatically selected when the thermometer is turned ON, follow these steps.

- Make sure the thermometer is turned OFF.
- Press and hold the blue SITE SELECTION button and, without releasing it, then press and hold the gray POWER button. Release both buttons when the site selection icon appears with 3 flashing dots.
- Press the SITE SELECTION button until the desired temperature measurement site is selected.
- Wait 4 seconds for the site indicator dot to stop flashing.
- Press the POWER button to turn the thermometer OFF.
- The preferred temperature measurement site will be automatically selected when the thermometer is turned ON.

**NOTE:** With a preferred temperature measurement site selected, an alternate temperature measurement site can still be selected during each use by pressing the blue SITE SELECTION button when the site selection icon appears after the thermometer is turned ON.

g. To remove the preferred temperature measurement site setting, open the battery door, remove the battery and reinsert it.

**IMPORTANT:** To ensure accurate readings, you must always take a temperature reading at the site on the individual that corresponds to the site selected on the thermometer (oral, underarm, or rectal).

### 5. PRECISION INSIGHT™ FEATURE

Precision InSight helps you to better understand the meaning of your child's temperature with the color indicated on the display. The screen displays Green if no fever, Yellow if fever, and Red if high fever. With Precision InSight™, the screen color is customized based upon the measurement site as shown on the chart to the right.

### 6. TURNING THE UNIT OFF

Turn the thermometer off by pressing the On/Off button. If you forget to turn it off, the thermometer will automatically turn off after approximately 30 seconds.

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